



## VIPM vip-file installation, (install OSA-driver in the LabVIEW instr.lib folder)

1. Try double clicking the vip file, if VIPM does not open:  
[download](#) and install VIPM, (free version).
2. If VIPM opens and quickly shuts down:
  - a. Open VIPM as administrator - either right click the VIPM exe-file or through Windows Start below.

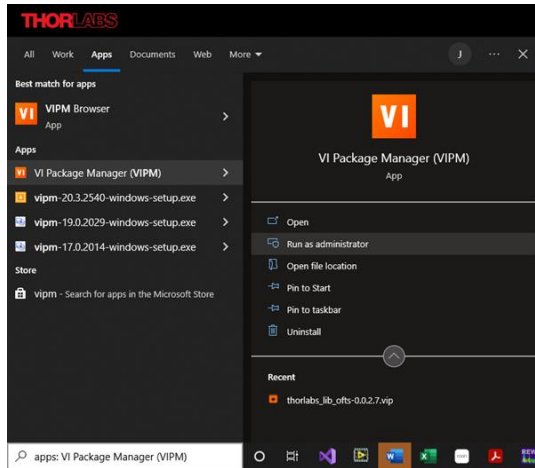


- b. Activate the menu and open the vip file, (File/Open Package File(s)..)

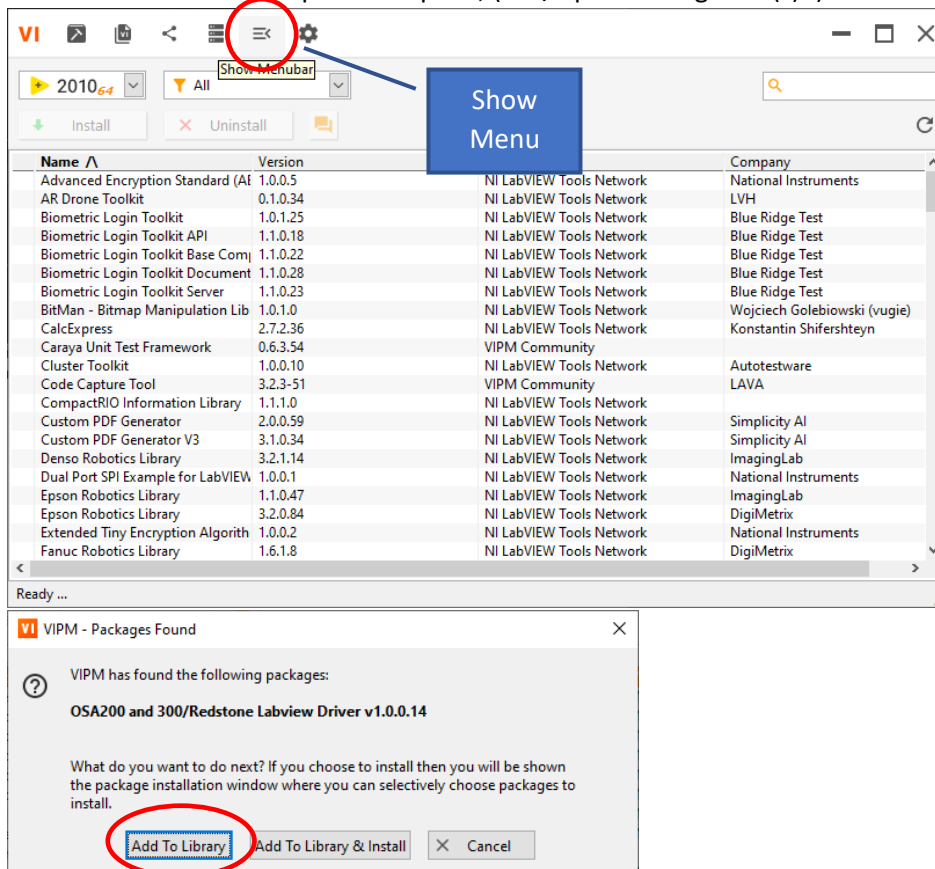


Continue next page..

3. Make sure the intended LabVIEW version is active/checked.

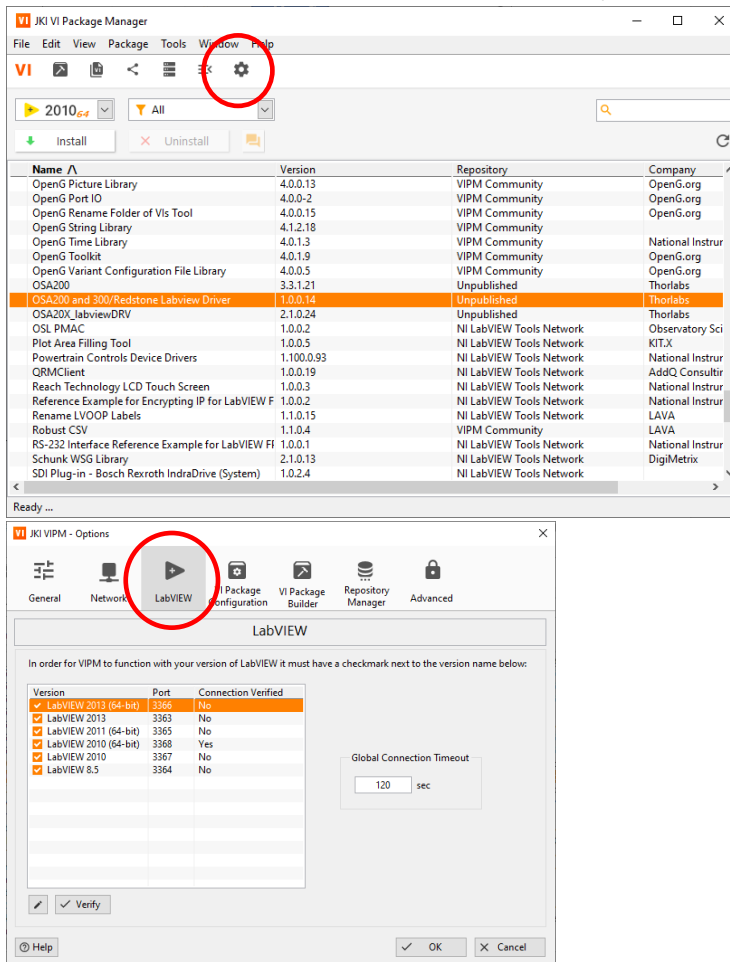


4. Double click the OSA-driver.

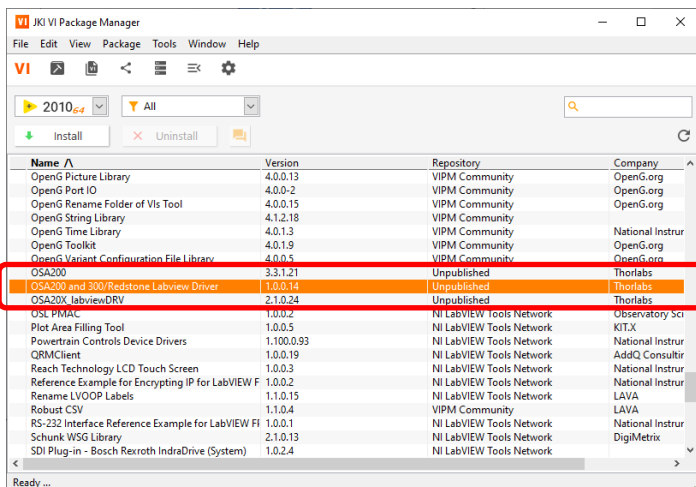

Continue next page..

5. Make sure no conflicting drivers are installed, (the combined OSA200 and 300 driver is not compatible with older OSA200 only drivers).

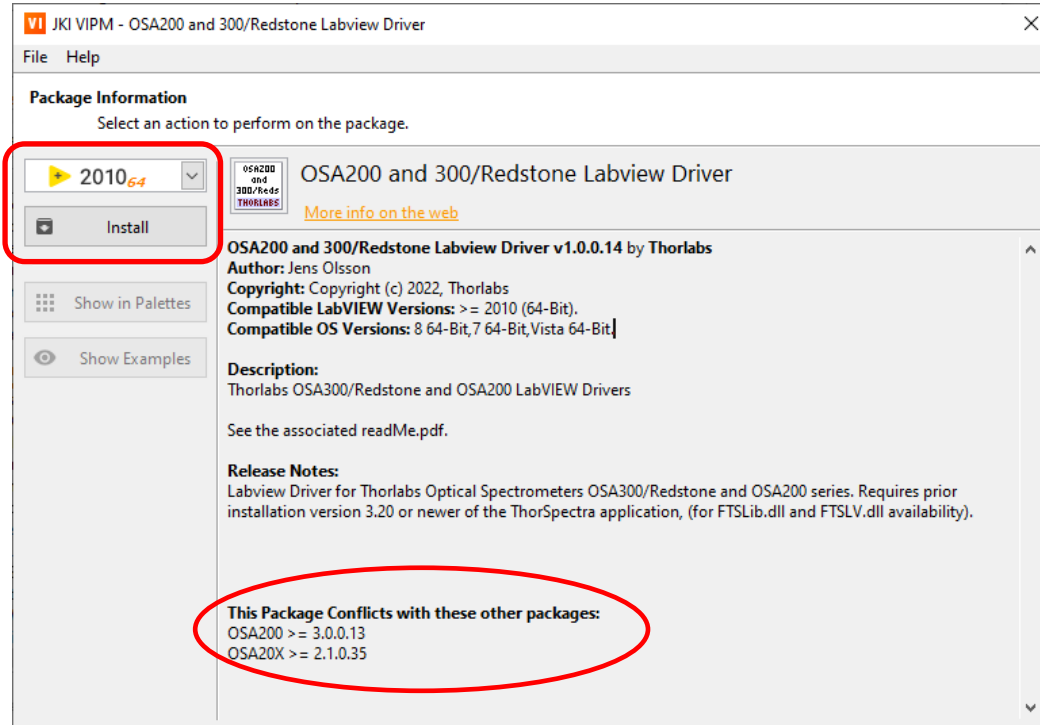

6. Select LabVIEW-version and install.